
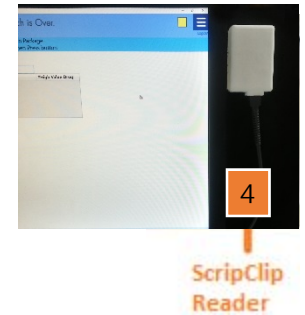
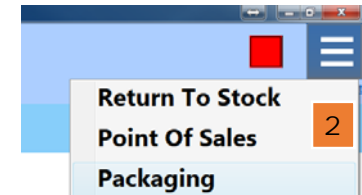


# Set Up

## Packaging Station:

- 1 Start the scripClip app.
- 2 Select *Packaging* from the pull-down menu
- 3 Select the color for the Packaging station by clicking on the color square of the scripClip application.
- 4 Position the Barcode scanner and the scripClip Reader as needed.

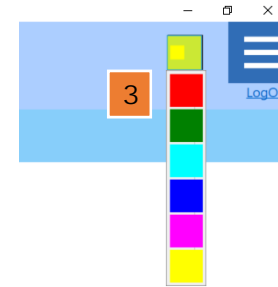

## Point of Sale Station:

- 1 Start the scripClip app.
- 2 Select *Point of Sale* from the pull-down menu
- 3 Select the color for the station by clicking on the color square of the scripClip application.
- 4 Position the scripClip Reader.

## Return Prescription to Stock:

- 1 Start the scripClip app.
- 2 Select *Return to Stock* from the pull-down menu
- 3 Select the color for the station by clicking on the color square of the scripClip application.
- 4 Position the scripClip Reader.

# Packaging the Prescription

## Hanging Bag scripClips:

- 1 Using the barcode scanner, scan one by one, up to 10 Rx numbers either from the vial or from the prescription paperwork for a customer.  
→ Scanned Rx numbers appear on the screen.
- 2 Place the vials and prescription paperwork in an unused Hanging Bag scripClip. Snap the hanging bag closed.
- 3 Hold the Hanging Bag scripClip near the scripClip Reader and press the orange button on scripClip.  
→ The scripClip will flash the station color once.
- 4 → scripClip application will display, 'Prescriptions clipped successfully.'



Place the filled hanging bag scripClip on the will-call rack.



# Packaging the Prescription

## Clip-on scripClips :

- 1 Using the barcode scanner, scan one by one, up to 10 Rx numbers either from the vial or from the prescription paperwork.  
→ Scanned Rx numbers appear on the screen.
- 2 Place the vials and prescription paperwork in the prescription bag.
- 3 Take an unused Clip-on scripClip **3a**, hold near the scripClip Reader **3b** and press the orange button on the scripClip.  
→ The scripClip will flash the station color once.  
→ scripClip application will display, 'Prescriptions clipped successfully.'
- 4 → scripClip application will display, 'Prescriptions clipped successfully.'
- 5 Attach the scripClip to a closed prescription bag. Place the prescription bag in the will-call area.

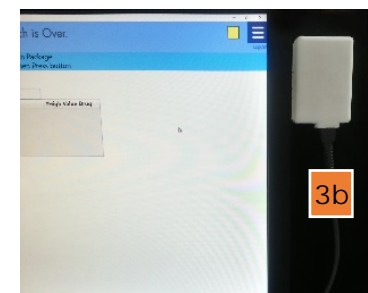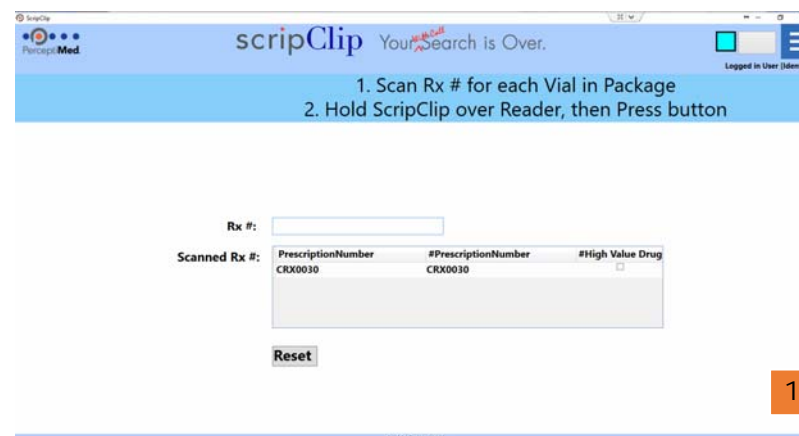

# Point of Sale: Locate a Prescription

## Locate a prescription:

- 1 Search for the Patient using name or date of birth (or both).
- 2 Press "Locate" to activate the scripClip's flashing light or "Locate with Tone" to activate the scripClip's flashing light and beep.
  - scripClip application shows "Located" and scripClip will flash for 30 seconds. **2a**
  - Retrieve the package with the flashing light that matches the station color.
  - Press the button on the scripClip to stop the flashing light.
  - Hold the Hanging Bag scripClip **2c** near the scripClip reader, and push the orange scripClip button. Or remove the Clip-on scripClip **2d** from the prescription bag, hold the scripClip near the scripClip Reader, and push the scripClip button
  - scripClip will flash the station color once to indicate transaction completed.
  - scripClip application will display 'Package Retrieved.' **2e**
  - Place empty scripClip in storage bin for use on next prescription package.



**2c**



**2d**

**1**

First Name	Last Name	Date of Birth	Phone	Zip	Pkg(s)
John	Miller	2/16/1966	(999)-555-2222	88445	1
John	Miller	8/26/1961	(999)-555-2222	88445	2

**2a**

Rx IDs	Status	Locate	Locate with sound
CRX0030	Located in Store	Locate	Locate with sound
CRX003P	Located in Store	Locate	Locate with sound

**2b**

Rx IDs	Status	Locate	Locate with sound
CRX001R	Package Retrieved	Locate	Locate with sound

**2e**

# Returning Prescriptions to Stock

## Return to stock:

- 1 Select number of days for “Prescriptions older than Days” from pull down menu. Press “Search”.
- 2 Scroll to bottom of list and press "Locate Packages" or "Locate Packages with Tone"
- 3 scripClip application shows “Located in Store” for each prescription on list and the scripClips flash for 30 seconds. Retrieve all of the packages with the flashing light that matches the station color.  
→ Press the button on the scripClip to stop the flashing light.
- 4 Hold the Hanging Bag scripClip near the scripClip Reader, and push the scripClip button (or) Remove the Clip-on scripClip from the bag, hold the scripClip near the scripClip Reader, and push the scripClip button.  
→ scripClip will flash the station color once to indicate transaction completed.
- 5 → scripClip application will display ‘Package Retrieved.’  
→ Place the empty scripClips in the storage bin for use on next prescription package.  
→ Repeat for each packaged to be returned to stock



Select “Packaging” from the pull-down menu in the upper right corner to return to the Point-of-Sale page.



# scripClip Tamper Alert (Admin Mode)

## Red Alert Message Issued:

- 1 Removing a scripClip from a package without using patient search or the Return To Stock screens causes scripClip to issue an alert. When the 'Alert' message appears, click on the *Alert* button, select the 'Locate' command to activate that scripClip.
- 2 After locating the scripClip, select the Action marked, 'Prescription OK' if the scripClip was accidentally removed or dislodged from the package. If the associated package indicates tampering or is missing, select the Action marked, 'Tampered/Missing Prescription'

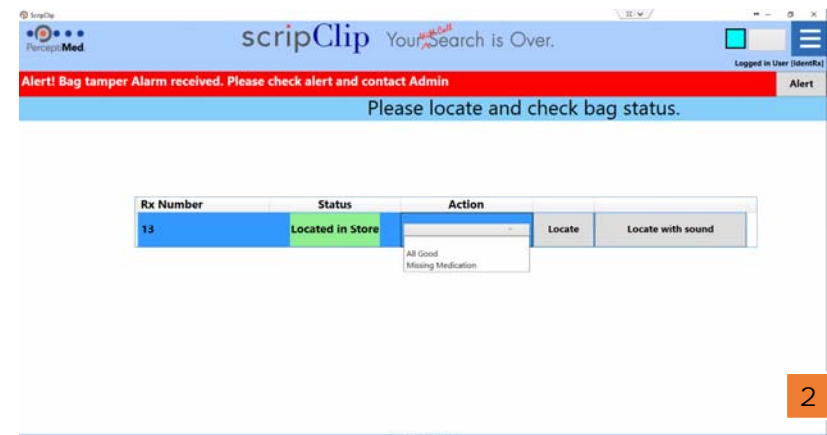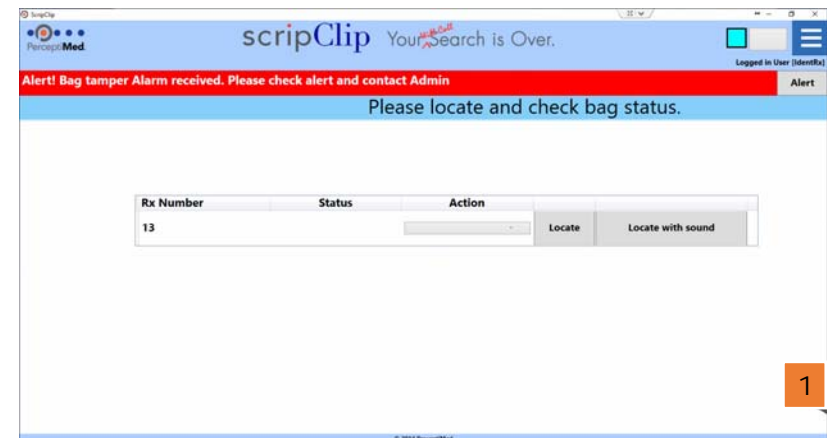

# Troubleshooting

Message: **Try again**

Cause: Package not located in store

**Solution:** Check Rx numbers entered. Press “Locate” or “Locate with Tone” again.

Message: **Wrong package retrieved.**

Cause: scripClip scanned does not match Rx numbers.

Solution: Put scripClip back on package and return to shelf. Press “Locate” again. Retrieve correct package.

Message: **Unable to read Rx numbers.**

Cause: scripClip reader did not receive all Rx numbers from scripClip

**Solution:** Hold scripClip closer to the scripClip Reader and press the orange button on scripClip again.

Event: **User is unable to find flashing scripClip in store**

Cause: scripClip stops flashing after 30 seconds.

**Solution:** Press “Locate Prescriptions” again or press “Locate Prescriptions with Tone”. User does not need to re-enter Rx numbers.

Message: **Unable to communicate with the central service**

Cause: Package not located in store

**Solution:** Call 650 941-7000 ext 5 or send email to

support@perceptiMed.com



# Appendix 1: Parts in the scripClip Kit

1 Hanging-Bag scripClip



2 Clip-on scripClip



3 scripClip Reader



4 scripClip Central Computers (x2)



5 scripClip Antennas (x4)





# Information



## Battery Replacement:

Do not attempt to replace the battery yourself. Perceptimed will replace battery and verify device as directed by service contract.

## Use Restriction:

Do not relocate system and components without a Perceptimed authorized installer. Components are restricted to indoor pharmacy and non-residential application.

## Disposal Information:

When this product reaches its end of life, return the device(s) to Perceptimed or dispose as directed by local laws and regulations. The collection and recycling of electronic waste conserves natural resources and protects human health and the environment.

# Compliance Information



## FCC Compliance Statement:

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is not recommended. If this equipment does cause harmful interference or is interfered with by other services contact Perceptimed or its representative.

Any modifications made to this device and accessories that are not approved by Perceptimed, Inc. may void the authority granted to the user by the FCC to operate this equipment.

## Applicable to the Following Devices:

Hanging-Bag scripClip, Models SCLP-HNGBAG-G1V1 and SCLP-HNGBAG-G1V2, FCC ID: 2AJ8D-HNGBAG-G1  
Clip-on scripClip, Model SCLP-CLPON-G1V1, FCC ID: 2AJ8D-CLPON-G1  
scripClip Reader, Model SCLP-READER-G1V1, FCC ID: 2AJ8D-READER-G1  
WAN Basestation, Model SCLP-WANTXRX1-G1V1, FCC ID: 2AJ8D-WANTXRX1-G1

# Compliance Information



## Disclosure of Exposure to Radiofrequency and Electromagnetic Fields:

This equipment complies with the FCC radiation exposure limits set forth per 47 CFR §2.1091 for an uncontrolled environment (ie: both occupational and unrestricted exposure to the public).

The operating power of the scripClip and its associated Reader is deliberately limited. No precautions are required for their use or operation.

The wide area broadcast station and its antennas must be installed and operated such that at least 20cm separation distance is maintained between these components and the operator or the general public. This is normally achieved by overhead ceiling mount of at least 8 feet from the level of the floor. These antennas may not be collocated with antennas for WiFi or other applications.